

Delivering F&P Optiflow™ + Tracheostomy High Flow Therapy



Introducing F&P Airvo™ 2, delivering respiratory support and airway hydration.

The Airvo 2 is an integrated solution for delivering Optiflow tracheostomy high flow therapy. Fisher & Paykel Healthcare's leading humidification technology allows the Airvo 2 to comfortably provide high flows of humidified air (with supplementary oxygen if required) to spontaneously breathing patients, through the unique Optiflow+ tracheostomy interface.

Optiflow+ tracheostomy interface

- F&P Evaqua™ technology, which reduces the formation of mobile condensate near the patient
- Options for both direct-connect and tracheostomy collar interfaces

Supplementary oxygen (when required)

- FiO₂ 21% to 100%
- Inbuilt ultrasonic oxygen analyzer requires no calibration, service or replacement

Designed for simple setup, use and cleaning

- Helpful onscreen animations assist with setup and troubleshooting

Adjustable flow settings

- Integrated flow generator delivers a wide flow range (10-60 L/min)
- No compressor or wall air supply required

F&P AirSpiral™ heated breathing tube

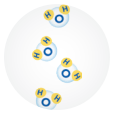
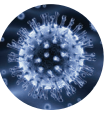


- No external probes, cables or adaptors required

Water vapor and other particles

The F&P Airvo 2 humidifier is designed to generate water vapor (molecules).

Nebulizers and other aerosol delivery systems are designed to generate liquid water particles. Unlike liquid water particles, water vapor is not able to act as a vector for transport of pathogens.

Sizes of particles and pathogens

| Pathogens | | | |
|---|---|---|--|
| Water vapor | Virus | Bacteria | Nebulized water particle |
|  |  |  |  |
| < 0.001 Microns | 0.02 – 0.3 Microns | 0.3 – 10 Microns | 1 – 10 Microns |

Product codes (*some accessories may not be available in all countries)

| | |
|-----------------------|---|
| PT101xx | Airvo 2 humidifier |
| 900PT561 | Airvo 2 AirSpiral tube and chamber kit (10-pack) |
| *900PT562, 563 | Airvo 2 AirSpiral tube and chamber kit with nebulizer adapter (10-pack) |
| OPT970 | Optiflow+ Tracheostomy interface (20-pack) |
| OPT971 | Sputum guard spare part (20-pack) |
| OPT980 | Optiflow+ Mask adapter interface (20-pack) |

Find out more at www.fphcare.com/airvo

F&P, AirSpiral, Airvo, Evaqua, and Optiflow are trademarks of Fisher & Paykel Healthcare Limited.

For patent information, see www.fphcare.com/ip

For further information about Optiflow+ tracheostomy interface including warnings and precautions, please read the user instructions.

613261 REV B © April 2021 Fisher & Paykel Healthcare Limited

Fisher & Paykel
HEALTHCARE